




Quest Full Shell Air Walker

Single Patient Use Only

Rx Only

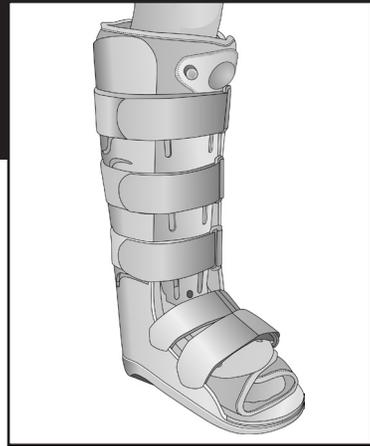

IMPORTANT INFORMATION

Please read all instructions before use. Correct application is essential for proper functioning of product.

INDICATIONS


Moderate to severe ankle sprains and strains; Stable foot and/or ankle fractures; Syndesmosis sprain; Post ORIF

PRECAUTIONS


This product is to be fitted by a physician (or properly licensed practitioner) who is familiar with the purpose for which they are responsible. The physician or practitioner is responsible for providing wearing instructions and precautions to other healthcare practitioners or healthcare providers involved in the patient's care and to the patient. Consult your physician (or properly licensed practitioner) immediately if you experience sensation changes, unusual reactions, swelling or increased pain while using this product. Follow the instructions of your physician (or properly licensed practitioner) for length and duration of use.

INSTRUCTIONS FOR USE


1. Unfasten strapping, remove front panel and open wrap.
2. Position foot and lower leg into wrap with the heel positioned in the back of the walker. Wrap liner around the foot then the leg and secure with hook fasteners.

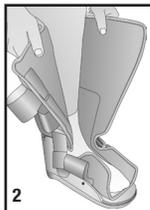

- 3a. Replace front panel.
- 3b. Secure strapping beginning with the foot and moving up the leg.

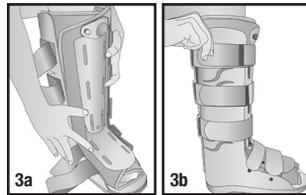

- NOTE:** Readjust as needed for a secure and comfortable fit.
4. To inflate liner, turn the inflate/deflate dial clockwise. Do not overtighten. Next, pump inflator bulb until desired compression is achieved.
 5. To deflate liner, turn the inflate/deflate dial counter-clockwise to release air pressure. Turn inflate/deflate dial clockwise when desired compression is achieved.
 6. Check for adequate circulation. Readjust air amount as needed for a snug and comfortable fit.

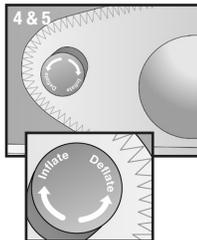

CARE INSTRUCTIONS

Wrap should be hand washed in cold water using mild detergent. Walker shell can be wiped down using mild detergent. Air dry.



Distributed by DeRoyal
200 DeBusk Lane
Powell, TN 37849 U.S.A.
888.938.7828
865.938.7828


Part #0-QuestWalker Revised 1/10
©2010 DeRoyal. All Rights Reserved.



Quest Full Shell Air Walker

Single Patient Use Only

Rx Only

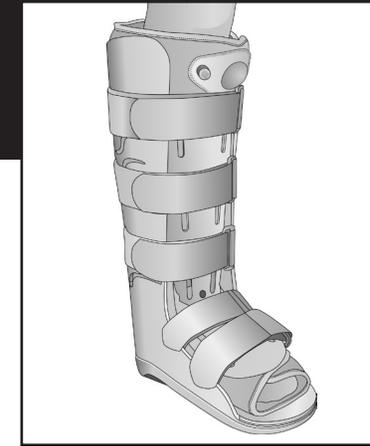

IMPORTANT INFORMATION

Please read all instructions before use. Correct application is essential for proper functioning of product.

INDICATIONS


Moderate to severe ankle sprains and strains; Stable foot and/or ankle fractures; Syndesmosis sprain; Post ORIF

PRECAUTIONS


This product is to be fitted by a physician (or properly licensed practitioner) who is familiar with the purpose for which they are responsible. The physician or practitioner is responsible for providing wearing instructions and precautions to other healthcare practitioners or healthcare providers involved in the patient's care and to the patient. Consult your physician (or properly licensed practitioner) immediately if you experience sensation changes, unusual reactions, swelling or increased pain while using this product. Follow the instructions of your physician (or properly licensed practitioner) for length and duration of use.

INSTRUCTIONS FOR USE


1. Unfasten strapping, remove front panel and open wrap.
2. Position foot and lower leg into wrap with the heel positioned in the back of the walker. Wrap liner around the foot then the leg and secure with hook fasteners.

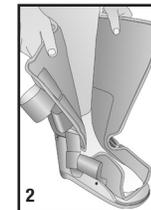

- 3a. Replace front panel.
- 3b. Secure strapping beginning with the foot and moving up the leg.

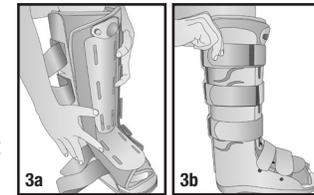

- NOTE:** Readjust as needed for a secure and comfortable fit.
4. To inflate liner, turn the inflate/deflate dial clockwise. Do not overtighten. Next, pump inflator bulb until desired compression is achieved.
 5. To deflate liner, turn the inflate/deflate dial counter-clockwise to release air pressure. Turn inflate/deflate dial clockwise when desired compression is achieved.
 6. Check for adequate circulation. Readjust air amount as needed for a snug and comfortable fit.

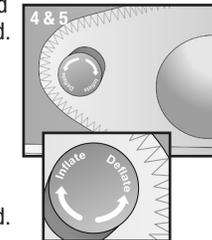

CARE INSTRUCTIONS

Wrap should be hand washed in cold water using mild detergent. Walker shell can be wiped down using mild detergent. Air dry.



Distributed by DeRoyal
200 DeBusk Lane
Powell, TN 37849 U.S.A.
888.938.7828
865.938.7828

Part #0-QuestWalker Revised 1/10
©2010 DeRoyal. All Rights Reserved.