



Orthoses are our passion

Elbow Brace Ulnaris

MR2250 · Elbow Brace Ulnaris

This elbow brace has been designed for ulnar nerve entrapments/cubital tunnel syndrome in the elbow. It is made from black ARX® ventilated material with Neptune™ material in the arm fold and features both a thermoplastic insert as well as a metal splint to prevent flexion in the elbow.


The thermoplastic insert and the metal splint can be used either separately or in combination for stability. The thermoplastic part should be molded externally by using hot water and then inserted into the pocket. Over the elbow joint there is a pocket on the inside with a memory foam pad that can be used to off-load a specific area by making a cut out in the foam pad.


The hook closure is made with low profile hook that can be attached anywhere over the elbow crease area on the orthosis. The design features a tongue of Neptune™ material to make the application easy.

Indications: Ulnar nerve entrapments, cubital tunnel syndrome, to limit flexion in the elbow joint for stroke patients or post-fracture treatments to limit elbow flexion.

Measurement chart

MR2250 · Elbow Brace Ulnaris

| Size | 1 (S/M) | 2 (L/XL) | 3 (XXL) |  |
|------|---------|----------|---------|---|
| cm | 20–25 | 25–31 | 31–38 | |

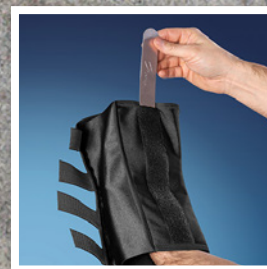

ARX® is a registered trademark by Mediroyal Nordic AB.

© Mediroyal Nordic AB 2021